

TempoCem®


Temporary Cements



All TempoCem temporary cements flow easily, for precise placement of cement and easy placement of the temporary.



Excess material can be easily removed with an explorer – especially important when cleaning the tooth for cementation of the permanent crown or bridge. The light color of the material makes it clearly distinguishable from teeth.

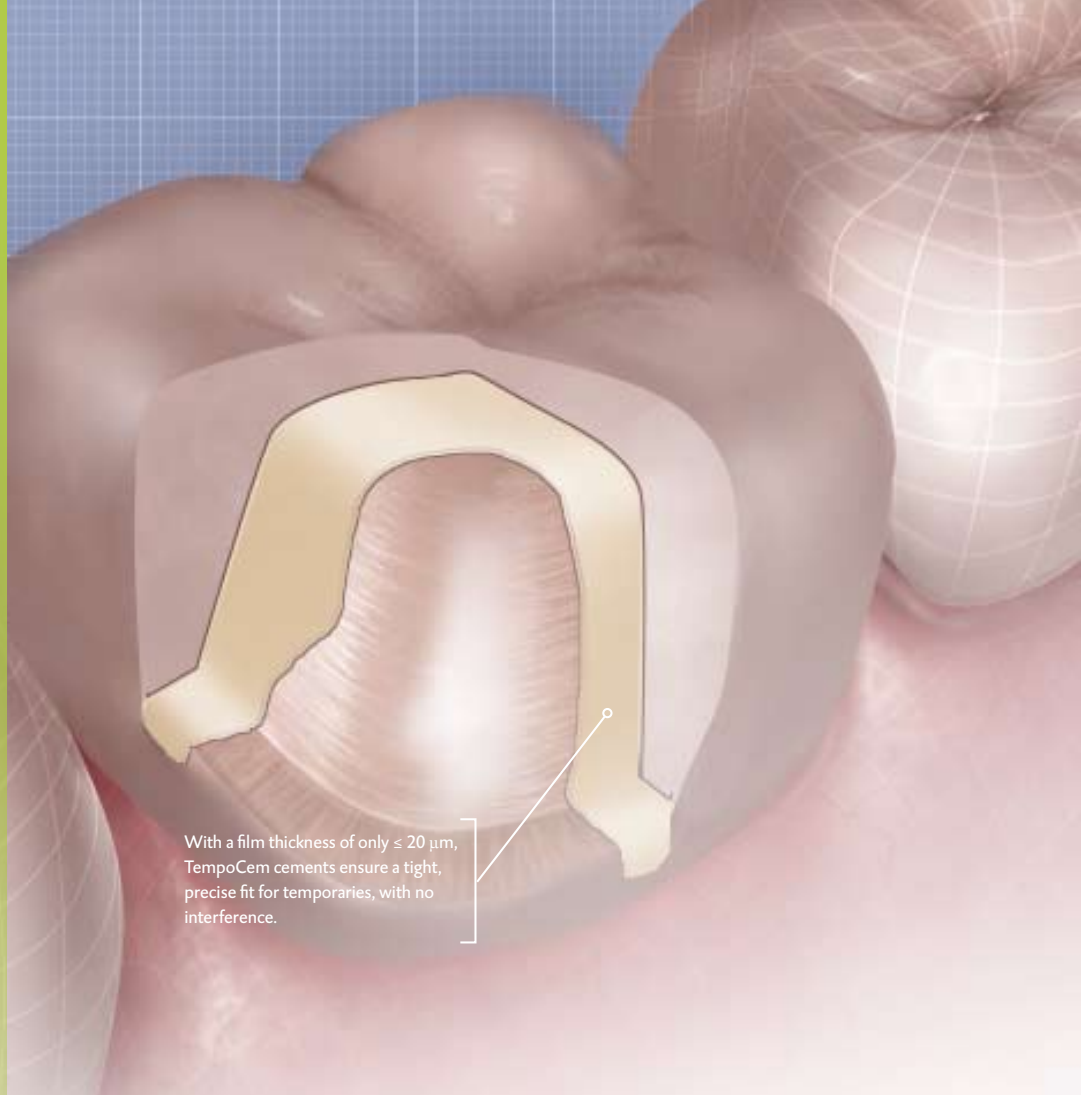

With a film thickness of only $\leq 20 \mu\text{m}$, TempoCem cements ensure a tight, precise fit for temporaries, with no interference.

TempoCem ensures an **accurate fit** for every temporary.

Three customized formulas provide reliable adhesion for all types of provisionals.

The zinc oxide/eugenol formula of original TempoCem temporary cement works as an analgesic for patients with extra sensitivity. Non-eugenol TempoCemNE is ideal for cases where the final restoration will be cemented with a resin or resin-reinforced glass ionomer. TempoCem Soft has the same chemistry as TempoCem, but sets to a softer end hardness for easier removal of tight-



fitting provisionals. All three are available in convenient, no-mess 25 ml Automix cartridges; TempoCem and TempoCemNE are also available in easy-to-handle Smartmix® syringes.

